

Anti-human α -Synuclein pAb 63

Catalogue number	Package size
0103000901	100 μ l
0103000903	1 ml


For research use only.

Product:	polyclonal antibody
Host:	rabbit
Immunogen:	recombinant human α -Synuclein
Purity:	undiluted serum
Storage:	at -20°C (repeated thawing and freezing should be avoided)
Applications:	ELISA, Western Blot (see figure below)
Dilution:	1:2.000-10.000 in ELISA, 1:1.000 in Western Blot
Reactivity:	human α -Synuclein

Figure

Recombinant human α -Synuclein was separated by SDS-PAGE under non-reducing (1) and reducing (2) conditions and blotted on nitrocellulose by Western blotting. NC membrane was incubated by pAb 63 1:10.000 diluted.

As detection system a horseradish peroxidase conjugated goat anti-rabbit IgG antibody and H_2O_2 /TMB were used. pAb 63 reacted with both α -Synuclein bands but with reduced protein much stronger.

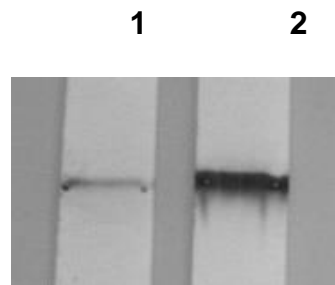

Roboscreen GmbH
Hohmannstraße 7
04129 Leipzig • Germany

Phone: +49 341 989 734 0
Fax: +49 341 989 734 199
info@roboscreen.com
www.roboscreen.com

CEO
Dr. Ingolf Lachmann

Ust-IdNr. DE 813064314
County Court Leipzig
HRB 17503
Leipzig