

*Anti-human TAU total pAb 64*

Catalogue number	Package size
0103001001	100 µl
0103001003	1 ml

*For research use only.*

---

<b>Product:</b>	polyclonal antibody
<b>Host:</b>	rabbit
<b>Immunogen:</b>	recombinant human TAU 441
<b>Purity:</b>	undiluted serum
<b>Storage:</b>	at -20 °C (repeated thawing and freezing should be avoided)
<b>Applications:</b>	ELISA, Western Blot (see figure below)
<b>Dilution:</b>	1:2.000-10.000 in ELISA, 1:2.000 in Western Blot
<b>Reactivity:</b>	human TAU


---

**Figure**

Recombinant human TAU was separated by SDS-PAGE under reducing (2) conditions and blotted on nitrocellulose by Western blotting. NC membrane was incubated by pAb 64 1:2.000 diluted.

As detection system a horseradish peroxidase conjugated goat anti-rabbit IgG antibody and H<sub>2</sub>O<sub>2</sub>/TMB were used. pAb 64 reacted with TAU protein.

**MW (kDa) 1**

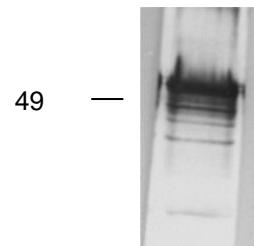

Roboscreen GmbH  
Hohmannstraße 7  
04129 Leipzig • Germany

CEO  
Dr. Ingolf Lachmann

Phone: +49 341 989 734 0  
Fax: +49 341 989 734 199  
[info@roboscreen.com](mailto:info@roboscreen.com)  
[www.roboscreen.com](http://www.roboscreen.com)

Ust-IdNr. DE 813064314  
County Court Leipzig  
HRB 17503  
Leipzig