

Anti-human tissue transglutaminase (E. coli) pAb R08

Catalogue number	Package size
0103000201	100 µl
0103000203	1 ml


For research use only.

Product:	polyclonal antibody
Host:	rabbit
Immunogen:	recombinant human tissue transglutaminase expressed in <i>E. coli</i>
Purity:	undiluted serum
Storage:	at -20 °C (repeated thawing and freezing should be avoided)
Applications:	Western Blot (see figure below), Immunofluorescence
Dilution:	1:1000 in Western Blot
Reactivity:	human tissue transglutaminase

Figure

Recombinant human tissue transglutaminase was separated by SDS-PAGE under reducing (2) conditions and blotted on nitrocellulose by Western blotting. NC membrane was incubated by pAb R08 1:1.000 diluted. As detection system a horseradish peroxidase conjugated goat anti-rabbit IgG antibody and H₂O₂/TMB were used. pAb R08 reacted with tissue transglutaminase protein.

MW (kDa) 1

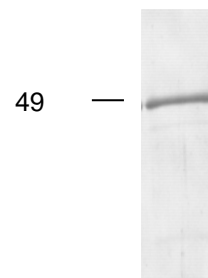

Roboscreen GmbH
 Hohmannstraße 7
 04129 Leipzig • Germany

CEO
 Dr. Ingolf Lachmann

Phone: +49 341 989 734 0
 Fax: +49 341 989 734 199
info@roboscreen.com
www.roboscreen.com

Ust-IdNr. DE 813064314
 County Court Leipzig
 HRB 17503
 Leipzig