

Anti-human tissue transglutaminase (insect cells) pAB R24

Catalogue number

847-0103000301
847-0103000303

Package size

100 µl
1 ml

For research use only.


Product:	polyclonal antibody
Host:	rabbit
Immunogen:	recombinant human tissue transglutaminase expressed in insect cells
Purity:	undiluted serum
Storage:	at -20 °C (repeated thawing and freezing should be avoided)
Applications:	Western Blot (see figure below), Immunofluorescence
Dilution:	1:1000 in Western Blot
Reactivity:	human tissue transglutaminase

Figure

Recombinant human tissue transglutaminase was separated by SDS-PAGE under reducing (2) conditions and blotted on nitrocellulose by Western blotting. NC membrane was incubated by pAb R24 1:1.000 diluted.

As detection system a horseradish peroxidase conjugated goat anti-rabbit IgG antibody and H₂O₂/TMB were used.

pAb R24 reacted with tissue transglutaminase protein.

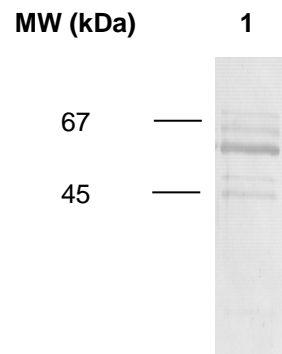

Roboscreen GmbH
Hohmannstraße 7
04129 Leipzig • Germany

CEO
Dr. Ingolf Lachmann

Phone: +49 341 989 734 0
Fax: +49 341 989 734 199
info@roboscreen.com
www.roboscreen.com

Ust-IdNr. DE 813064314
County Court Leipzig
HRB 17503
Leipzig