

Anti-human β -amyloid mab RS-07

Catalogue number

 847-0102008501
 847-0102008503

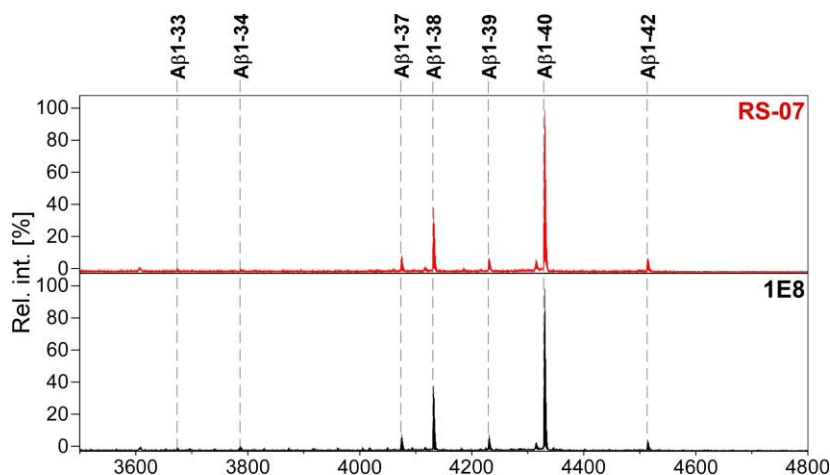
Package size

 100 μ g
 1 mg

For research use only.

Product:	monoclonal antibody
Clone:	RS-07
Host:	mouse
Isotype:	IgG1 _k
Immunogen:	Synthetic peptide VKMDAEFRHDSGY
Specificity:	β -amyloid N-terminal epitope, 1-38, 1-40, 1-42
Purification:	purified by affinity chromatography
Purity:	> 95% by SDS-PAGE
Solution:	PBS pH 7.4 without additives
Storage:	at -20 °C (repeated thawing and freezing should be avoided)
Applications:	ELISA, Western Blot, Immunohistochemistry
Dilution:	2 ng/ml in ELISA, 1 μ g/ml in Western Blot
Reactivity:	β -amyloid N-terminal

IP-MS from pooled human CSF using RS-07 coupled on M270 expoy beads (Thermo Fisher, USA) compared to commercially available antibody 1E8 (Nanotools, Germany). Analysis performed and results transferred kindly by Prof. Dr. Olaf Jahn and Dr. Hans Klafki, University Goettingen.


 Roboscreen GmbH
 Hohmannstraße 7
 04129 Leipzig • Germany

 CEO
 Dr. Ingolf Lachmann

 Phone: +49 341 989 734 0
 Fax: +49 341 989 734 199
info@roboscreen.com
www.roboscreen.com

 Ust-IdNr. DE 813064314
 County Court Leipzig
 HRB 17503
 Leipzig